


SINGLE-USE CARTRIDGE: ANOTHER TRUPOINT ADVANTAGE

An important feature of the IRMA TRUpoint system, is the single-use cartridges which are available in a variety of analyte configurations. Each cartridge contains miniaturized electrochemical measurement sensors and an onboard calibrant. With a wide operating test range, IRMA TRUpoint can be used in colder ORs. Most cartridges do not require refrigerated storage and are ready to use anywhere, anytime.

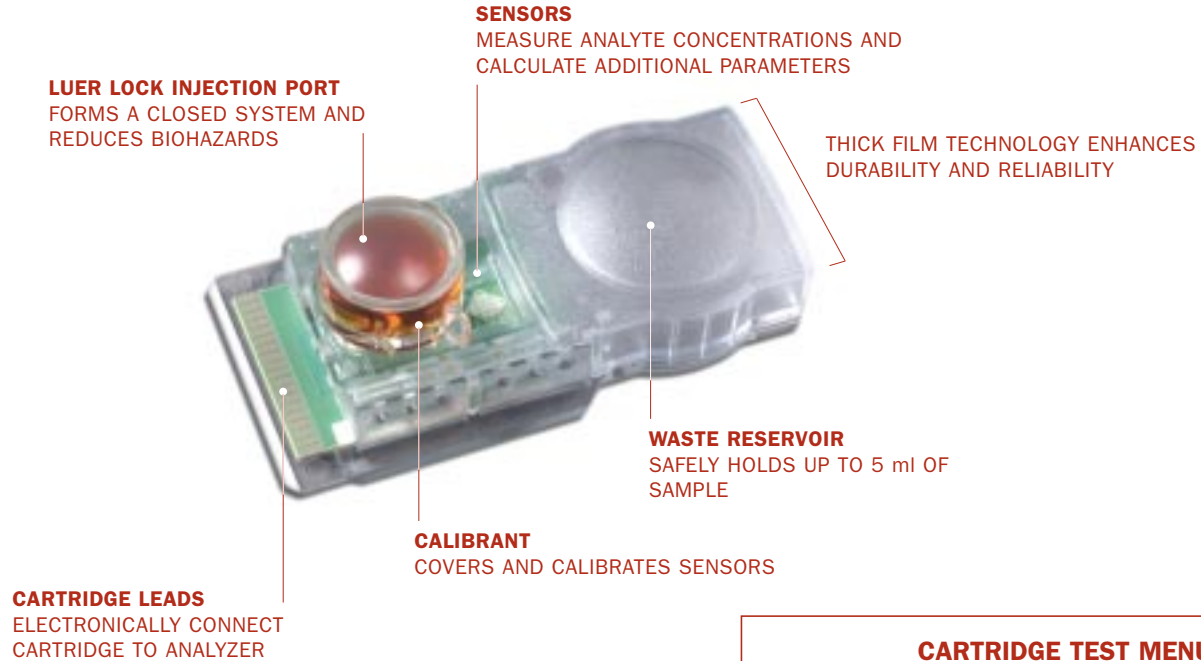

		CC	BG	H3	H4	GL	CR
Measured Values	pH	•	•				
	pCO ₂	•	•				
	pO ₂	•	•				
	Na ⁺	•		•	•	•	
	K ⁺	•		•	•	•	
	Cl ⁻				•	•	
	iCa	•		•			
	Glu					•	
	BUN				•		
	Hct	•		•	•		
Creatinine						•	
Calculated Values	HCO ₃ ⁻	•	•				
	TCO ₂	•	•				
	BEecf	•	•				
	BEb	•	•				
	O ₂ SAT	•	•				
	iCa(N)	•					
	tHb	•		•	•		

ITC


For more than 35 years, ITC has been at the Point of Care. By focusing exclusively on bedside testing longer than anyone else, we've always made a difference wherever you need it. And with service and support that's second to none, we're by your side, so you can stay by theirs.

To learn more about IRMA TRUpoint and how ITC can help you improve patient care at the bedside, contact us at 1.800.631.5945 or visit www.itcmed.com/irma.shtml.

©2004 ITC. All rights reserved. IRMA TRUpoint is a trademark and IRMA is a registered trademark of ITC.



8 Olsen Avenue • Edison, NJ 08820 USA
tel: 732.548.5700 • fax: 732.248.1928
ITC Europe • Strada Rivoltana • 20090 Rodano (MI) ITALY
tel: +39.02.9532.0196 • fax: +39.02.9532.0276
www.itcmed.com
a subsidiary of Thoratec Corporation



The True Point of Care.




Accurate results where and when you need them. IRMA TRUpoint sets the standard for testing blood gas, electrolytes and other chemistries right at your patient's side—the true point of care. This powerful, self-contained system can easily be transported from OR to CCU to ICU, or wherever you need it for immediate results and more timely diagnostic interventions.

Bedside testing that's accurate,

ACCURATE

Obtain immediate results right at the bedside for improved efficiencies, diagnostic capabilities and procedural outcomes.

- Precise results within two minutes
- Lab-quality accuracy
- Comprehensive menu
- Custom correlation for site-specific analyte reference ranges

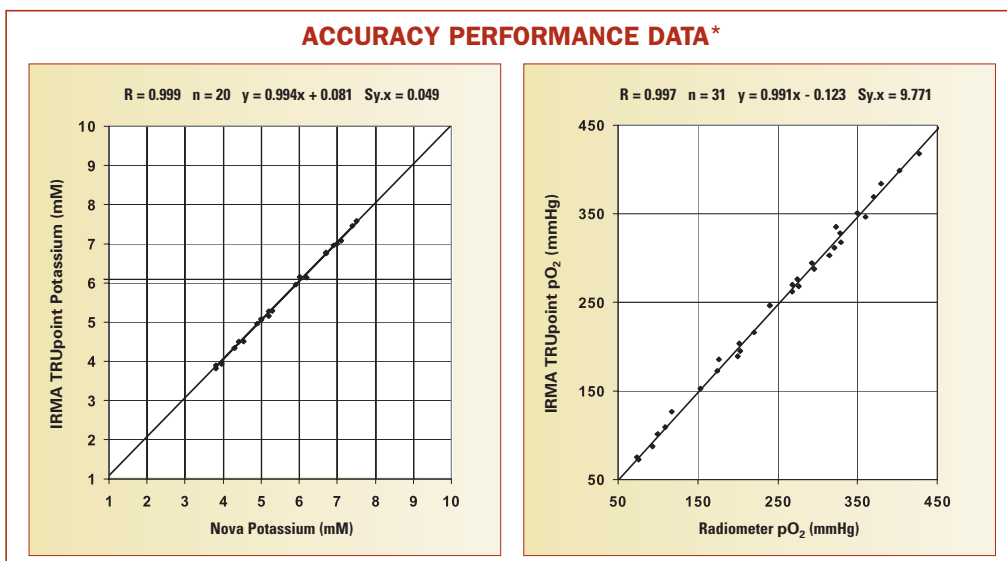

* Verification studies at multiple sites

ONBOARD THERMAL PRINTER



OPTIONAL
BAR CODE SCANNER



INTERACTIVE TOUCH SCREEN

SINGLE-USE CARTRIDGE



reliable and self-contained.

RELIABLE

Trust, integrity and performance with simple, intuitive operation and advanced quality features.

- Single-use cartridge for maintenance-free convenience
- Self-contained calibrant in each cartridge
- Real-time, automatic electronic quality control for enhanced safety and improved compliance
- Exceptional service and support

SELF-CONTAINED

Discover complete portability with lightweight, easy-grip, battery-operated, all-in-one design.

- Interactive touch screen
- Onboard thermal printer documents results and test notes
- Complete data management from patient information to lot traceability
- Bar code capabilities for on-the-spot accuracy
- Comprehensive connectivity, including modem, serial port or Ethernet



IRMA
TRUPOINT[™]
Blood Analysis System