

Retrieval Technology

LoopMaster-Sochman Snare



- **First Choice Safety Device For Special Foreign Body Retrieval**
- **High Performance Guaranteed With NITINOL Loop**
- **Highest Product Safety With Atraumatic Tip**


Especially designed by the Team of
Dr. Sochman & Dr. Peregrin in Prague
For most difficult retrieval manoeuvres



DELIVERING LESS INVASIVE FUTURE TECHNOLOGY


LoopMaster-Sochman Snare

1. Atraumatic distal guide wire tip
2. Nitinol Snareloop
3. Complement guidewire
4. Radiopaque platinum inner marker
5. Highly flexible shaft
6. Supportive sleeve
7. Luer connector
8. / 9. Extra "BIGG" torque device to ease manipulation

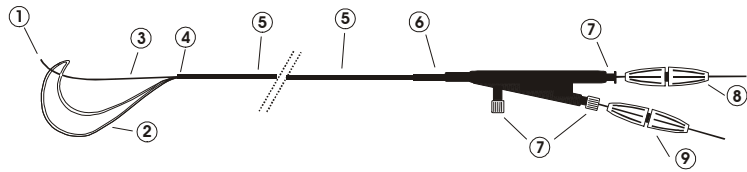

* The Loopmaster-Sochman Snare has been specifically designed to retrieve foreign bodies without an open end. It makes it your first choice safety device for the most difficult procedures.

Retrieval Technique Of Non Open End Foreign Bodies

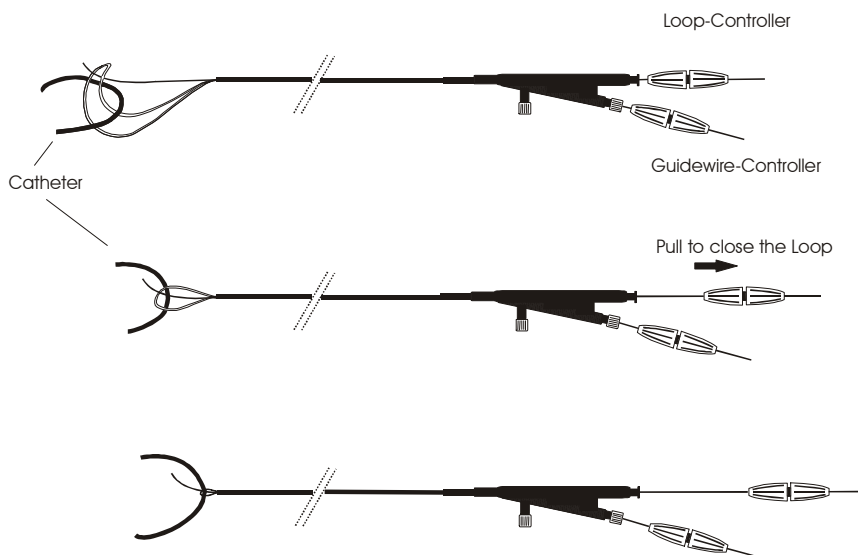

1. Open the Loop and manipulate around the foreign body

2. Closing the loop will allow to catch the foreign body between snare and guidewire

3. Safe retrieval of the foreign body with closed loop

Order Information Loopmaster-Sochman Snare:

<u>Product Code</u>	<u>Loop-Diameter (mm)</u>	<u>Length (cm)</u> <u>Loop/GW</u>	<u>Introducer</u> <u>Sheath (cm)</u>	<u>Introducer</u> <u>Sheath (F)</u>
ASL-25	25 mm	130/150	110	8

High Quality Products Made In Germany!

Your local Distributor:



Andramed GmbH
 Silvanerstraße 8
 D - 72766 Reutlingen
 Germany
www.andramed.com
 Tel.: +49-7127-972290
 Fax: +49-7127-972291

CE 0483
 Rev. 01/07